

BugStopper™

Sterile venting for culture vessels

BugStopper is a novel device that gives you a simple and reliable way of venting your culture vessels. This reusable, sterile closure, is made of biosafe silicone rubber and incorporates a hydrophobic glass microfiber filter that provides an aseptic vent.

Available in two sizes, BugStopper will fit most widely used culture flasks. BugStopper allows air to enter and exit the flask while excluding bacteria and viruses. This eliminates contamination of both the contents as well as the outside environment.


The self-healing silicone wall of the BugStopper can be pierced with needles for inserting a sample port or for gas infusion. This is not recommended for the BugStopper 10 because of the wall thickness.

The BugStopper is easy to use; you simply push it into the neck of the flask and autoclave the entire vessel. It is also economical because it can be autoclaved again and again without loss of performance. GE Healthcare has certified 50 cycles.

BugStopper 10 has an enhanced feature, a protective shroud covering the lip of the flask.

Materials of construction

BugStoppers are manufactured from biosafe silicone rubber. The vent is a composite of hydrophobic ultrafine glass microfiber, reinforced with polyester monofilament laminates. A stainless steel reinforcement ring surrounds the vent for added support.

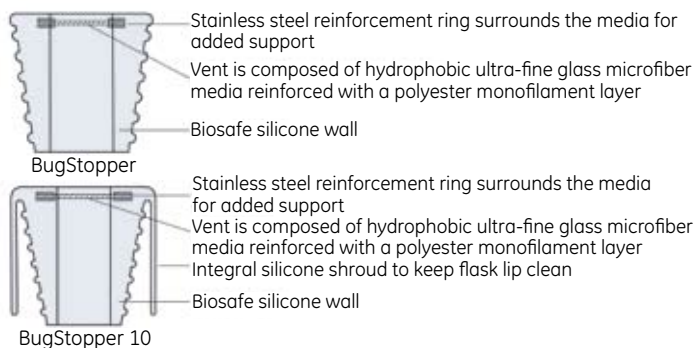

Easy ordering procedure

The Whatman™ BugStopper is available through our extensive network of national distributors. You can purchase any of these products in the same manner as your other lab supplies and receive the superior service to which you are accustomed.

Ordering information

Catalog number	Description	Quantity
6713-3010	BugStopper	10
6713-3100	BugStopper	100
6713-6010	BugStopper 10	10
6713-6050	BugStopper 10	50

For more information, visit www.gelifsciences.com/whatman

GE Healthcare Bio-Sciences Corp.
800 Centennial Avenue
Piscataway, NJ 08855-1327 USA



GE, imagination at work and GE monogram are trademarks of General Electric Company. Whatman and BugStopper are trademarks of GE Healthcare companies. All third party trademarks are the property of their respective owners. © 2005-2009 General Electric Company – All rights reserved. Previously published Nov. 2005.

All goods and services are sold subject to the terms and conditions of sale of the company within GE Healthcare which supplies them. A copy of these terms and conditions is available on request. Contact your local GE Healthcare Representative for the most current information.