




GD/X

Filters samples with high particulates with less hand force

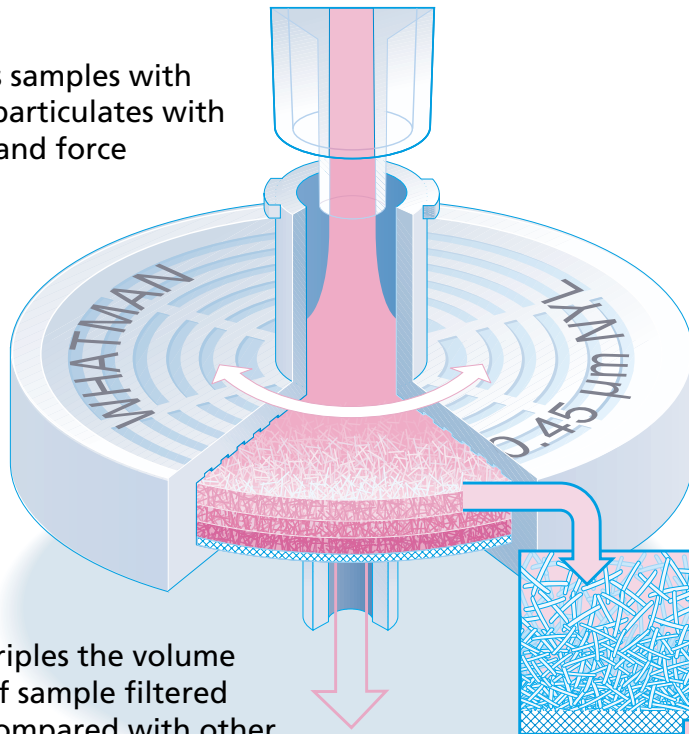

Triples the volume of sample filtered compared with other conventional filters

Extends filter life and reduces costly exchanges

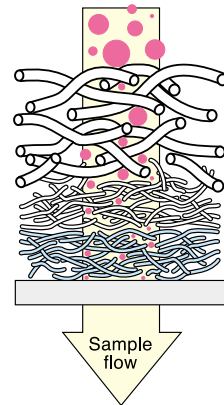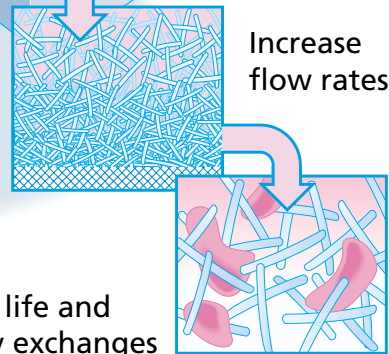

GD/X syringe filters contain four filtration layers which substantially reduce blockage and increase volume throughput



Layers 1 & 2
Glass microfiber prefilter
GMF 150 (graded density)
filters from 10 μm to 1 μm



Layer 3
Glass microfiber prefilter
GF/F filters down to 0.7 μm



Layer 4
Membrane filter
Nylon, PVDF, PTFE, etc.
Filters down to 0.45 μm or 0.2 μm